


BOUGIE

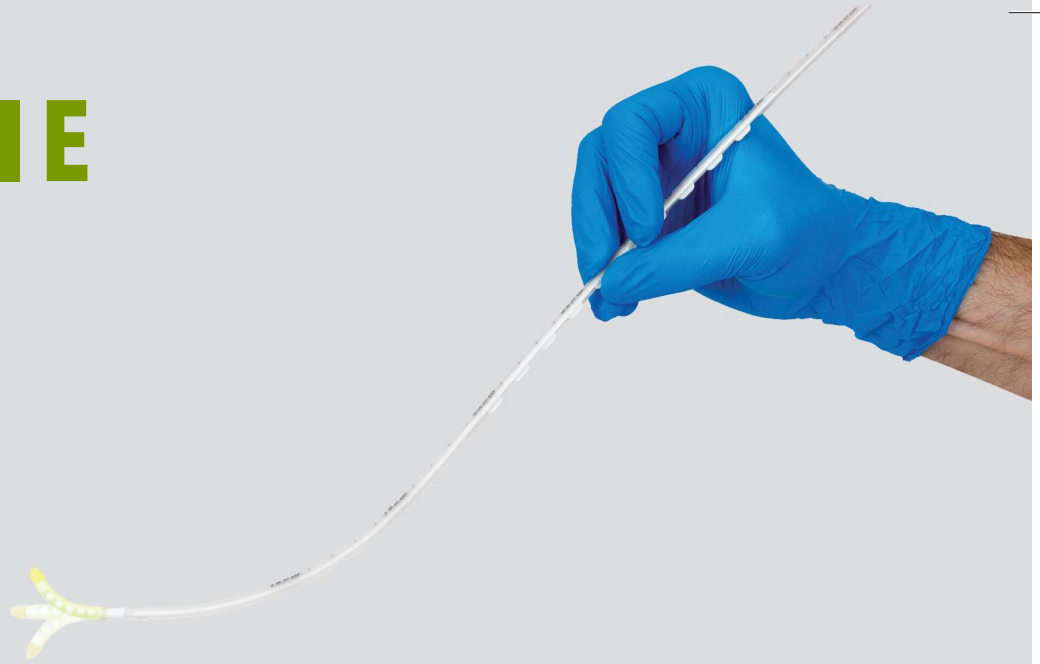

FLEXIBLE TIP BOUGIE

- Flexible, steerable tip makes it easier to navigate difficult airways
- Innovative slider is grasped between the thumb and forefinger to move the tip anteriorly and posteriorly
- Blunt silicone tip reduces the risk of tracheal trauma
- Pre-curved design to provide anterior flexion, use the control mechanism to retro-flex and straighten to to move past the vocal cords
- Slider tabs Slider tabs prevent the Bougie from passing through the vocal cords
- Bright phosphorous coating on tip makes it highly visible, especially when used with Intubrite dual ultra-violet and LED lights
- Sterile single-use device
- Tabs on slider are easy to grip and manoeuver particularly when used with gloves
- Use in conjunction with video laryngoscope for optimal results
- Graduation marks from 10cm to 50cm (65cm overall length)
- Use with ETT's with 7mm inner diameter or larger

NAVIGATING THE AIRWAYS IS EASIER WHEN YOU HAVE CONTROL

Bougie introducers are an invaluable tool when dealing with challenging airways, but sometimes they're just not enough. The flexible, steerable Tip Bougie makes it easier to navigate difficult airways. It's blunt silicone tip reduces the risk of tracheal trauma and the bougie's pre-curved design enables the control mechanism to retroflex and straighten to move past the vocal cords, thus making it safe to use.

THE FLEXIBLE TIP BOUGIE WILL HELP YOU FIND THE WAY

Simply grasp the flexible Tip Bougie between your thumb and your forefinger, and then use your thumb to move the integrated slider forward and backward motion to flex the tip anteriorly and posteriorly to navigate the device between the cords and into the larynx.

